

AIVIA SOLO

Standalone
AED docking
station



AIVIA

Protection and autonomous monitoring

The Aivia SOLO docking station is fully energy self-sufficient and uses a patented technological innovation. AIVIA SOLO monitors and protects the AED while maintaining the medical device within the manufacturers recommended operable temperature range.

To meet the need for autonomy of the AED docking station in outdoor public places, the AIVIA SOLO innovates by facilitating deployment with a simple ground installation without any electrical connection.

The use of a defibrillator in the first five minutes following cardiac arrest increases the chances of survival.

The installation of standalone Aivia SOLO terminals that are visible and identifiable in the urban environment makes it easy to locate the nearest defibrillator, act very quickly and save lives.



Energy autonomy

No connection
to the electricity grid
Recyclable materials
Energy stored



Remote monitoring

AED presence and operability
AED environment
Door opening and/or closing
Geotracking



Deployment

Easy to install
Simple ground layout
Consult us for your
deployment projects





Technical specifications

Monocrystalline Solar Panel	48W
Heater	40W
PCM HC8 phase change materials	Heat capacity 208 Wh
Insulating thermal envelope	1,17 K.m ² /W
Dimension	L 602cm x H 2400cm
Weight	60kg
Connection	IoT Sigfox, NB-IoT, 4G-LTE



Monitor all of your AEDs

The AIVIANet platform allows management of all the AEDs contained in surveillance cabinets.

The platform gives a view on the AED indicator and all the information relating to its environment (opening, removing the AED, temperatures outside the ranges, expiration of consumables, etc.).



Intuitive

Viewing alerts
Geotracking of AEDs
Optimized user experience



Security

Securing information
Materiovigilance
Access rights management




Configurable

Configuring alerts and recipients
Web service for interfacing with external software



Monitoring

Operational management
Automatic alerts by email or sms

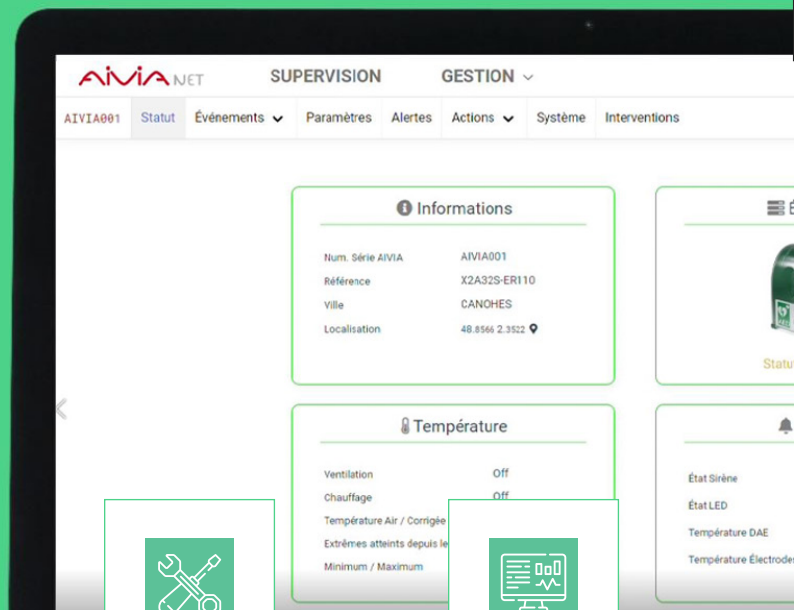



Hayes First Aid

Supplies | Training | Equipment

www.hayesfirstaid.com

t. 069 66006 | 01 5079900

e. sales@hayesfirstaid.com

w. www.hayesfirstaid.com