

iPAD SP1 Defibrillator & Heated Outdoor Cabinet

A complete outdoor-ready defibrillator package — the device, the heated locked cabinet, the rescue kit and the signage. Everything you need to put a defibrillator on an outside wall and keep it working.

WHAT YOU ARE ACTUALLY BUYING

Ten years of not having to think about it.

Anyone can sell you a defibrillator in a box. The box is the easy part. The hard part is the ten years afterwards — the pads that expire quietly, the battery that runs out in year five, the person who looked after it leaving in year three and taking everything they knew with them.

That is the part we do. We fit it out ready to use, we hold the dates, we tell you before anything expires rather than after, and there is one named person here who knows your site and answers the phone.

You are not buying a device from us. You are handing the whole ten years to somebody whose job it is to remember.



A HAYES FIRST AID INSTALL OF THIS CABINET
— PERMANENT-LIGHT VERSION SHOWN



Adult & child

ONE SET OF PADS AND A SWITCH
— NO CHILD PADS TO BUY



10 years

WARRANTY WHEN YOU REGISTER
IT
7 YEARS AS STANDARD



Built for outdoors

HEATED, LOCKED, IP66 CABINET
RATED FOR IRISH WEATHER



The defibrillator

iPAD SP1 Semi-Automatic — supplied with battery, dual-use adult/child pads and carry case, fitted and tested before dispatch

CU-SP1-S



The cabinet

Heated, IP66, code-locked polycarbonate clamshell with a 10-year enclosure guarantee — the standard version, our best seller

DFS4K442LLE-Y

Product images are the manufacturer's / distributor's own. The install photograph is a Hayes First Aid installation of the same cabinet type, not your site.

01

KEY FEATURES & SPECIFICATIONS

The iPad SP1 is a semi-automatic defibrillator. It analyses the casualty's heart rhythm, decides whether a shock is needed, charges itself, and then tells the rescuer to press the orange shock button. It will not allow a shock on a rhythm that does not need one.



One switch for adult or child
The pads are dual-use — no second set to buy, find or fit under pressure.



It gets louder when it needs to
Ambient noise detection raises the voice prompts to 90 dB, so instructions carry over traffic or a crowd.



CPR detection and a metronome
It paces the rescuer rather than just counting at them.



Robust one-piece design
Nothing to unclip, break or lose. Survives a 1.2 m drop to any edge or corner.



Status shown on the face
Pad life, battery level and overall readiness — readable without opening anything.



Event data on an SD card
Up to 5 events, 3 hours each, downloaded with the device still on the wall. No sending it away.

DEVICE SPECIFICATION — MODEL CU-SP1	
AED model	CU Medical Systems iPad SP1 Semi-Automatic — SKU CU-SP1-S
Operating mode	Semi-automatic (rescuer presses the shock button)
Defibrillation waveform	e-cube biphasic, truncated exponential type
Output energy	150 J at 50 Ω load for adults · 50 J at 50 Ω load for children
Time from CPR to shock	Less than 6 seconds from the completion of CPR to shock delivery
Charging time	11 seconds typical, on a new battery pack
Patient impedance range	25 Ω to 175 Ω
Guidance	Voice prompts and illuminated visual instructions, with ambient noise detection to 90 dB
Adult / child pads	SMART dual-use adult and child pads — SKU CUA1007S. Switchable at the device; no separate paediatric pads required. Typical 24-month shelf life — the date printed on the packaging governs.
Battery	Lithium-ion, non-rechargeable, 12V 4.2Ah — SKU CUSA1103BB. 5 years from installation, minimum 200 shocks.
Dimensions & weight	260 × 256 × 69.5 mm (W × L × H) · 2.4 kg including battery and pads

Sealing	IP55 (IEC 60529)
Drop protection	Withstands a 1.2 m drop to any edge, corner or surface
Operating conditions	0°C to 43°C, humidity 5–95% non-condensing, 0–15,000 ft
Event memory	5 events, up to 3 hours per event, SD card download with free software
Warranty	7 years as standard · free upgrade to 10 years when you register the device

What the energy figures mean in practice

Resuscitation guidance is built around minimising the pause in chest compressions — the time the casualty has no circulation. The figure that matters on this device is **under 6 seconds from the end of CPR to shock delivery**, and the dual-use pads remove the one step most likely to cost time on a child: finding and fitting a different set of pads.

Energy and timing figures are the manufacturer's stated nominal specifications for model CU-SP1. Clinical first-shock efficacy data for this model is not held on file at Hayes First Aid and is therefore not quoted here.

02 WHAT IS IN THE BUNDLE

Bundle SKU CU-SP1-SPBO — 4 items

€1,700.00 · 0% VAT on the defibrillator and cabinet



The defibrillator

IPAD SP1 · CU-SP1-S



The cabinet

HEATED, LOCKED, IP66



BLS rescue kit

FREE — DIRECT ORDERS



AED wall sign

FREE — DIRECT ORDERS

List price shown. Your quotation carries the price actually offered to you and takes precedence over this sheet.

ITEM	DETAIL	LIST
IPAD SP1 Semi-Automatic CU-SP1-S	Supplied rescue-ready — battery and dual-use adult/child pads fitted and tested before it leaves us.	€1,095.00
Heated polycarbonate outdoor cabinet, locked DFS4K442LLE-Y	Double-skinned clamshell in impact-resistant polycarbonate — UV-stabilised and fire-retardant to UL94 V0. IP66 , dust-tight and water-resistant. 465 × 435 × 240 mm (H × W × D), 10 kg. Thermostatic 50W heater. Internal hook and mounting kit included. This is the standard version — a permanent-light version is available separately if you want the window illuminated. 10-year enclosure guarantee. Fits all AED brands. Requires a mains supply — see section 03.	€620.00
Green Pouch BLS CPR Rescue Kit 5001049G	Gloves, face shield, shears, razor and wipes — the things you need in the first minute and will not want to go looking for. Supplied free with every defibrillator on a direct order.	€13.50

AED Heart Restarter sign, rigid
S3218R

So people can find it. **Supplied free with every defibrillator on a direct order.**

€8.00

Not included: electrical installation and the 240V supply to the cabinet (section 03) · wall fixings suited to your particular wall · a spare set of pads or a spare battery, if you want one on the shelf.

03 INSTALLING IT — WHAT YOU NEED TO ARRANGE



Registered Electrician Required — Before Delivery

This cabinet requires a **240V mains power supply**. Electrical installation must be carried out by a registered electrician in accordance with Irish wiring regulations (ETCI). Hayes First Aid supplies and delivers the cabinet only — electrical installation is not included in the price.

- 1 Engage a registered electrician to prepare a 240V mains supply point at your intended mounting location
- 2 We deliver the cabinet — you mount and connect it with your electrician
- 3 Once powered, the heating system is thermostatically controlled — no ongoing adjustment needed

Not sure what's involved? Call us on **069 66006** — we will send you the manufacturer's installation instructions to give your electrician.

Choosing where it goes

- **Somewhere findable, and reachable at night.** A defibrillator nobody can locate in the dark is not available. Lighting matters as much as height.
- **Think about who will actually fetch it.** The nearer the front of the site, the better.
- **Mounting height** should suit the people who will use it, and the cabinet needs to open fully without fouling a door, downpipe or sign.

Wall type decides the fixings. The cabinet ships with a mounting kit; masonry, timber frame, cladding and insulated render all behave differently, and that is a question for whoever mounts it.

The cabinet is code-locked — so here is how the code is looked after

A code-locked cabinet only works if the code exists in more than one place. That is the whole point of recording it properly.

We record the cabinet, not just the defibrillator. At onboarding we take the cabinet's **make, model, serial number and access code** — for every locked cabinet. That is part of how we look after it: our engineer has to be able to open it at the annual visit, and the information then exists somewhere other than in one person's memory.


It is never published. The code sits on your record for our engineers and for you. Your Guardian record will only ever say that the cabinet is locked and to call **112 or 999** — a code is never shown to a passer-by.

Register with the National Ambulance Service too (section 05). Their form asks whether the cabinet is locked and, if it is, asks for the access code — so that a 112/999 call handler can pass it to whoever is running to get it.

And tell more than one person on site. Whoever looks after it should never be the only one who knows — which is why we hold a primary guardian plus two backups as a minimum.

04 YOUR NEXT TEN YEARS, IN FULL

Two things expire on a known clock and the device tells you when it is unhappy. Here is every event between now and 2036, what it costs, and which of us is doing it.

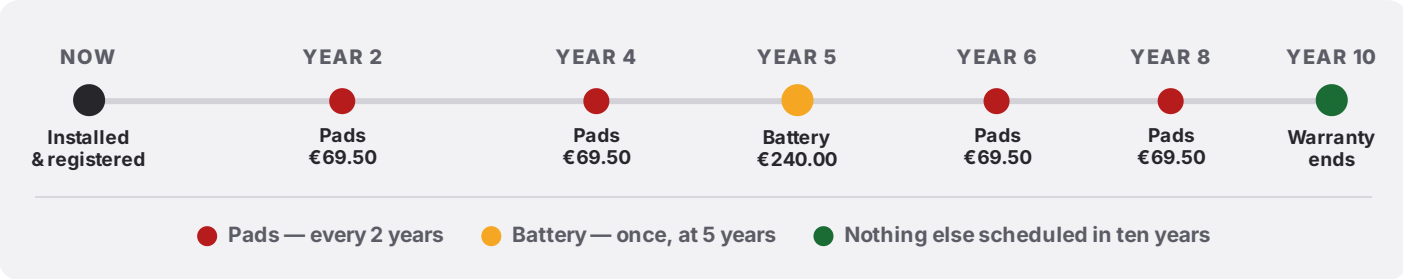

YEAR	WHAT HAPPENS	WHAT YOU DO	COST
Now	It arrives with the battery and pads already fitted and tested, in its cabinet, with the rescue kit and the sign. Your device goes on the Guardian register and your reminders are set.	Get the electrician in, mount it, register it, tell us who the guardians are	€1,700.00
1	On Guardian Managed, an engineer comes to you, inspects it, tests it, checks the dates, updates your record.	Let us in and show us where it is	included in Managed
2	Pads are due. We tell you at 90, 60 and 30 days, and we have them on the shelf.	Just say yes	€69.50
3	Annual visit. And if the person who looks after it has left, the record has not.	Nothing at all — that is rather the point	included in Managed

YEAR	WHAT HAPPENS	WHAT YOU DO	COST
4	Pads again. Same three reminders.	Just say yes	€69.50
5	Battery is due — five years from installation. We have been counting since the day it shipped.	Just say yes	€240.00
6	Pads.	Just say yes	€69.50
7	Annual visit. Reminders. The record keeps itself current.	Nothing at all — that is rather the point	included in Managed
8	Pads.	Just say yes	€69.50
9	Annual visit.	Nothing at all — that is rather the point	included in Managed
10	Warranty runs to here if you registered it. We will tell you what your options are before you have to ask.	Nothing at all — that is rather the point	—
Ten years	Equipment and consumables — the device, the cabinet, four sets of pads and one battery		€2,218.00

That is the **complete** list of things that will happen to the equipment — no annual renewal on the equipment, no hidden consumable, nothing you have to remember. Pads are costed on a 24-month cycle and the battery on its five-year life from installation; both carry 0% VAT. The dates printed on the pad packaging and the battery govern in practice — a pad set fitted late in its shelf life expires sooner. Prices are today's list prices and are not held for ten years.



You cannot order the wrong pads

This device takes one specific pad set and nothing else — SKU **CUA1007S**. Your model and serial number are on your record, so when a reminder goes out it is for the right part, and when you order we check it against your record before anything ships. **You never have to know the part number.**

WHAT THE DEVICE DOES ON ITS OWN, AND WHAT IT WANTS FROM YOU

Self-testing	The SP1 tests itself and shows the result on the status indicators on its face — readiness, battery level and pad life.
Your routine check	Look at the indicators. That is the whole job. Hayes First Aid recommends a documented check at least monthly, and we supply the checklist.
If an indicator changes	Ring us on 069 66006 with the model and serial number to hand. Do not open the device.
After it has been used	Ring us the same day. The pads are single-use and must be replaced, the battery may need replacing depending on how many shocks were delivered, and the event data on the SD card can be downloaded and converted.

Disposal

Used pads and spent batteries are hazardous waste — they do not go in general waste. Hayes First Aid operates a take-back service; civic amenity sites also accept them free of charge.

05

REGISTERING IT — THREE REGISTRATIONS, ALL YOURS TO DO

REGISTER WITH	WHY	WHO DOES IT
National Ambulance Service National AED Network	An unregistered defibrillator is invisible to a 112/999 dispatcher. Registering means a call handler can send someone to <i>your</i> device. It is also where the locked-cabinet access code is recorded.	You register directly — it's free and takes minutes. Start at www.hayesfirstaid.com/register-your-aed
CU Medical device warranty	Registering the device upgrades the warranty from 7 years to 10 years , free.	You register it. The details come in the box and we will point you at the right page.
AED Guardian Hayes First Aid	So the expiry dates, the serial numbers, who looks after it and every check done to it live somewhere other than one person's memory. See section 06.	You register it — and on Guardian Managed we sit down and build the record with you at onboarding, so it is a conversation rather than a form.

Paediatric use

Covered as standard. The SP1's SMART pads are dual-use — a switch on the device changes it to child mode, delivering 50 J instead of 150 J. **There are no separate paediatric pads to buy, store or remember.** Optional separate child pads exist but are not required.

Current resuscitation guidance is to use a defibrillator on a child rather than withhold one. If in doubt in an emergency, follow the device's prompts and the instructions of the 112/999 call handler.

06

AED GUARDIAN — CHOOSE YOUR LEVEL



AED
GUARDIAN
by Hayes First Aid

Guardian is the part of Hayes First Aid that remembers. **Every defibrillator we supply goes on Guardian at no charge** — the dates are held and the expiry reminders start before the box reaches you. From there you choose how much of the work you want to keep, and how much you want us to take off you.

Guardian

Included

FREE WITH EVERY AED WE SUPPLY
DIRECT

WHAT YOU GET

- Your device on the register — **model, serial, install date, pad and battery expiry dates**
- **Expiry reminders at 90, 60 and 30 days** — these can never be switched off
- Primary guardian plus two backups
- A dated record you can ask us for at any time
- Your record held and kept current by us
- QR stickers in the onboarding pack — you apply them
- No annual visit at this level
- No AED Inspection Record at this level

Guardian Pro

€5 /month

€60 A YEAR EX VAT · ONE AED

EVERYTHING IN GUARDIAN, PLUS

- **A named Hayes First Aid contact** for queries — a person, not a ticket number
- Your own dated check record, kept for you
- Additional AEDs at **€60 a year each**
- Still self-managed — no engineer visit
- No AED Inspection Record at this level
- No consumables discount at this level

For sites happy to do the checking themselves but who want somebody to ring.

MOST POPULAR — RECOMMENDED

Guardian Managed

€15 /month

€180 A YEAR EX VAT · ONE AED

EVERYTHING IN GUARDIAN PRO, PLUS

- **We come to you once a year.** A trained engineer inspects the device, verifies the pad and battery dates and works the full checklist on site
- **The AED Inspection Record** — every check, visit and replacement, dated and kept
- **10% off every consumable** for as long as you subscribe — €6.95 off each set of pads, €24.00 off each battery on this device
- Your engineer applies and tests the QR stickers at the first visit
- A named contact who **knows your site and its history**
- Additional AEDs at the same site at **€120 a year each**

Rolling subscription — billed annually in advance, no fixed term to sign, and it keeps going until you tell us to stop. Prices are ex VAT and are per defibrillator at one site. Ten or more devices, or several sites? That is priced for what you actually have — ask us. **You can change level at any renewal.**

Why Managed is the one we recommend for this install

Four reasons, and the first one is specific to an outdoor cabinet on the mains.



More can go wrong outdoors

The heater, the seal, the code, the hinges, the fixings — and none of it is visible from an office. Once a year somebody who knows what he is looking at stands in front of it.



Nobody inspects a cabinet from a screen

Any service implying otherwise is worth a hard question. Managed is the level where **the visit is real.**



The 10% pays for a chunk of itself

Ten years on this device is four sets of pads and one battery — **€51.80** back on equipment you were buying anyway.




The record stops being anyone's job

The person responsible today will move on. The dates, the history and the reminders **stay exactly where they are.**

A missed check doesn't just sit there

If a check is missed we remind you, then your backup guardians — and if the defibrillator still hasn't been checked, we pick up the phone. **Five steps before we ever let it rest, and the fifth one is a person.**

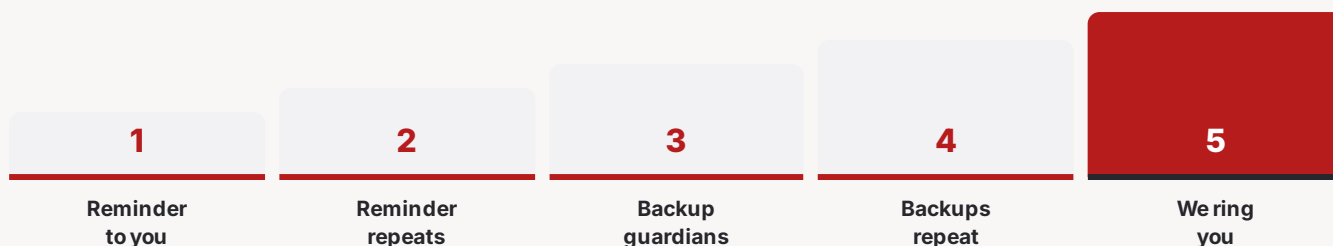

You set it up the way that suits you — and you can change it whenever



How often we visit

Once a year is the default and it is what most sites want. If you would rather we came twice a year, or every quarter, say so and we will price it. Your call, not ours.



Who the guardians are

Add people, remove people, change the primary. We hold a primary plus two backups as a minimum, because the one person who knows must never be the only person who knows. When someone leaves, tell us and the reminders follow the new name.



How often you check it yourself

We recommend **the interval the manufacturer states for your device** — we set that on your record and send you the checklist. You can tighten it if you want to; checking more often than we ask is never held against you.



How the reminders reach you

Mobile message, email, or both — set per person, so the site manager can have texts and head office can have email. **Two kinds of reminder, and only one is optional:** the *time to check your AED* nudge can be turned down or off if you find it a nuisance. The **pad and battery expiry reminders cannot be turned off** — those keep coming at 90, 60 and 30 days, whatever else you switch off.



It rolls

No fixed term to sign. It renews each year and continues until you tell us to stop.

What happens when you sign up to Managed

1. **We build the record with you.** Model, serial number, install date, pad and battery expiry dates, the cabinet's make, model, serial and access code, and the names of your primary guardian plus two backups.
2. **You get a confirmation** with what we hold and what is due when. Check it — if a date is wrong, that is the moment to say so.
3. **The reminders start.** 90, 60 and 30 days before anything expires, to your named people, by the channel each of them chose.
4. **We book your first annual visit** and work around your site. The engineer inspects, tests, checks the consumables, updates the record and leaves you with it.
5. **Your named contact is introduced** — a person, with a direct number.
6. **It renews once a year** until you tell us otherwise.

What Guardian is not, so there are no surprises

Guardian is **not insurance and not a compliance certificate**, and it will never tell you that you are compliant — for most of this there is no legal standard to be compliant with. What it does is say **what you have, what state it is in, and what is due**. The authority comes from the record being true.

It is also **not a substitute for registering with the National Ambulance Service** (section 05). That registration is what puts your defibrillator in front of a 112/999 dispatcher, and it is separate, free, and yours to do.

Works with all major AED brands — Philips, HeartSine, Zoll, Mediana, Mindray, LIFEPAK and more. Any AED. Any brand. Any supplier.

07

WHY BUY IT FROM US

You can buy this defibrillator from a website that will never speak to you again. Here is what is different about buying it here.



It arrives ready to save a life

Battery in, pads on, self-test passed, before it leaves us. You are not assembling a medical device on an outside wall in the rain.



We start counting the day it ships

Most people find out a pad expired when they open the cabinet. You will find out ninety days before.



A person, not a ticket number

One named contact who already has your history in front of them. You don't explain yourself from the beginning every time.



The record outlives the person

Whoever is responsible today will move on. The dates, the history and the checks stay exactly where they are.



We come to you

On Managed an engineer arrives once a year and does the whole checklist himself. Nothing is posted anywhere and the device never leaves your wall.



We tell you when it's finished

When this one reaches the end of its life we will say so early and tell you your choices, rather than letting you find out when the parts stop being made.

Trading since 2012 — and still here to answer for it

Hayes First Aid has been supplying, training and maintaining since 2012. We do all three, which is the reason the maintenance half of this document is as long as the specification half. A supplier who only sells has no reason to care what happens in year five.

Next working day delivery to all 32 counties. Order before 12pm Monday to Friday for same-day dispatch, subject to stock.

Rescue-ready

TESTED BEFORE
DISPATCH

0% VAT

ON THE AED AND
CABINET

32 counties

NEXT WORKING DAY

One person

WHO KNOWS YOUR
SITE

The 0% VAT position

Defibrillators, their electrode pads, their batteries and wall-mounted cabinets are **zero-rated for VAT in Ireland** — VATCA 2010, Schedule 2, paragraph 11(3)(ba). That is why the two main items in this bundle show 0%. Guardian subscriptions are a service and are quoted ex VAT.

Hayes First Aid

Castle Demesne, Newcastle West
Limerick, V42 N220, Ireland
Mon–Fri 09:00–17:00

T 069 66006 · 01 507 9900

E sales@hayesfirstaid.com

W www.hayesfirstaid.com

MacroMed Healthcare Limited

t/a Hayes First Aid
Registered in Ireland No. 510559
VAT IE9824728C

iPAD SP1 & Heated Outdoor Cabinet Bundle · sales sheet v1.9 · 3 September 2026 · bundle SKU CU-SP1-SPBO. Free-with-defibrillator items and the Guardian inclusion apply to direct orders; reseller and trade orders are quoted separately. Specifications are the manufacturer's stated nominal figures for model CU-SP1. Prices are list prices at the date of this sheet and exclude electrical installation. Guardian subscription prices are ex VAT. Where this sheet and your quotation differ, your quotation governs. In an emergency call 112 or 999. Prepared by Hayes First Aid.